

THE VERSATILE CHOICE



Birthright converts quickly and easily from a standard bed to a versatile delivery platform, allowing a wide choice of safe and comfortable positions to be adopted. The counterbalanced base provides optimum stability.

PRODUCT SPECIFICATIONS

Product Description	Order Code
BIRTHRIGHT with mattress, safety sides, all accessories* and accessory storage trolley with cover	170000EL
BIRTHRIGHT with mattress, safety sides, all accessories* and wall-mounted accessory storage rack	170001EL
BIRTHRIGHT with mattress and safety sides but no accessories	170010EL
BIRTHRIGHT with mattress but no safety sides or accessories	170011EL
Accessories	Order code
Leg rest (single)	170103
Strap for leg rest (single)	35104
High foot support, left	170105*
High foot support, right	170106*
Low foot support, left	170107*
Low foot support, right	170108*
Hand grip (single)	170109*
Fluid collection bowl	170110*
Step	170111*
IV drip rod	170112*
Leg rests with straps (pair)	170114*
Safety sides (pair)	170115*
Accessory Storage Trolley	170100
Accessory trolley cover	170113
Accessory storage rack (wall-mounted)	679.120

(* denotes accessories included as standard)

Mattress Re-Order Codes	Order Code
Mattress including cover - body section	170101
Mattress including cover - foot section	170102
Cover only - body section	170101/COVER
Cover only - foot section	170102/COVER

Technical Data	
Overall length	2210 mm (87")
Overall width	1000 mm (39.25")
In bed length normal	2050 mm (80.75")
In bed length retracted	1380 mm (54.25")
Height range (to top of mattress)	660 - 960 mm (26 - 37.75")
Safe Working Load	180 kg (396 lbs)
Product weight (bed only)	130 kg (286 lbs)
Backrest angle	0-65°
Trendelenburg Tilt angle	0-12°

Electrical Data	
Power in	3A max. @ 230V a.c. 50-60Hz 5.8A max. @ 120V a.c. 60Hz
Power out	24V d.c. max. 70VA
Electric shock protection	Class 1, Type B
Ingress protection	IP54
Electrical Safety Standards	Complies with EN60601-1

ArjoHuntleigh is a branch of Arjo Ltd Med. AB.
Only ArjoHuntleigh designed parts, which are designed specifically for the purpose, should be used on the equipment and products supplied by ArjoHuntleigh. As our policy is one of continuous development we reserve the right to modify designs and specifications without prior notice.
® and ™ are trademarks belonging to the ArjoHuntleigh group of companies.
© ArjoHuntleigh, 2009



Standard bed position



Semi-chair position with high foot supports and safety sides raised



Chair position with handgrips and step



Chair position with low foot supports and safety sides raised



Chair with high foot supports and handgrips



Lithotomy position with leg rests

ARJOHUNTLEIGH
GETINGE GROUP

www.ArjoHuntleigh.com

ArjoHuntleigh
310-312, Dallow Road Luton
Bedfordshire LU1 1TD
United Kingdom
Phone: +44 1582 745 800
Fax: +44 1582 745 866
E-mail: international@arjohuntleigh.com

GETINGE GROUP is a leading global provider of products and systems that contribute to quality enhancement and cost efficiency within healthcare and life sciences. We operate under the three brands of **ArjoHuntleigh**, **GETINGE** and **MAQUET**. **ArjoHuntleigh** focuses on patient mobility and wound management solutions. **GETINGE** provides solutions for infection control within healthcare and contamination prevention within life sciences. **MAQUET** specializes in solutions, therapies and products for surgical interventions and intensive care.

Printed in the Czech Republic • Point • 01.MB.14.1.GB-INT.1.AHG • December 2009

BIRTHRIGHT
DELIVERY BED RANGE



...with people in mind



THE PRACTICAL CHOICE



Convenient storage: wall-mounted accessory rack


Accessories and storage

The *Birtright* system includes a full range of attachable accessories, designed to optimize positioning for delivery or clinical procedures.

A mobile trolley keeps accessories neatly stored when not in use. Alternatively, a wall-mounted accessory rack is available for areas with limited space.

A SINGLE SOLUTION FOR LABOUR, DELIVERY AND RECOVERY

Birtright™ is an electrically operated delivery bed designed to provide support for mother, and efficient functionality for nursing staff, throughout all stages of the birthing process. Its flexibility and ease-of-use offers a practical, cost-effective solution for the modern, progressive labour suite.

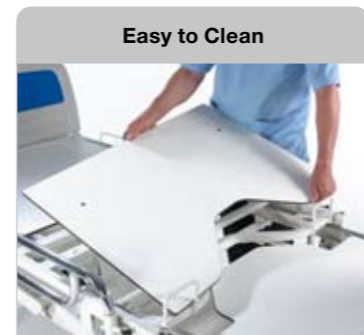


Removable mattress platform sections make the bed easy to clean.



Electrically operated height and backrest adjustments reduce handling risks.



Dual-sided levers lower the backrest instantly in the event of an emergency.

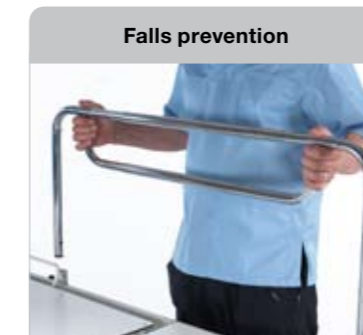

Adjustable safety sides, fitted as standard, are removable when not required.

Additional features and benefits:

- Manual Trendelenburg tilt to 12° for clinical management
- Retractable platform with 2-part mattress, optimizing access to the perineal area
- Removable platform sheets for easier cleaning and improved infection control
- Dual-side manual CPR release levers for backrest
- Central brake/steer castor system for optimum stability and manoeuvrability
- Battery back-up in the event of mains power failure
- Comprehensive range of easy-to-fit accessories